



Traveller testing with mariPOC®

Trusted diagnostics partner offering rapid and affordable COVID-19 solutions

The high-capacity mariPOC® COVID-19 antigen test is well suited for traveller testing at transport hubs such as airports, harbors, train stations and bus terminals. The test swiftly detects the infectious individuals enabling **safe travels**. The COVID-19 super-spreaders can be identified in 20 minutes. With mariPOC even the individuals with low amounts of virus but contagious are identified.

Key features

- 300 tests per day with 1 mariPOC analyser
- Market-leading sensitivity and specificity
- Detects different virus variants
- Automated result readout
- No lab training required to perform the test
- Only one person is needed for analysing the samples

Performance

Sensitivity (N)

Study 1:

PCR Ct ≤ 30: 100% (13/13)

Study 2*:

PCR Ct ≤ 28: 100% (31/31)

PCR Ct ≤ 30: 91.9% (34/37)

PCR Ct ≤ 34: 84.4% (38/45)

Specificity (N)

100% (201/201)

Reference test

PCR

*NOTE! Study 2 relied on nasopharyngeal swab samples in transport medium (UTM/VTM). Compared to using native swabs, transport media samples can affect adversely to the stability of samples and the performance of antigen tests. The results of Study 2 should be interpreted accordingly.

Order information:

Ref: 1204S | mariPOC® SARS-CoV-2

Tests/plate: 308

Testing capacity: 308 samples/24 h

Easy work flow



Nasopharyngeal swab sample



1 min sample handling



Fully automated analysis



Most positive results in 20 min*

* Negative results confirmed in 55 min.




Workplace testing solutions

Ensuring the health and safety of workers and customers is a top priority

The high-capacity mariPOC® COVID-19 test offers solutions tailored to your requirements, helping you keep your operations open in a safe and manageable way. The test swiftly detects the infectious individuals helping you maintain a safe work environment. mariPOC can be installed and operated at the company's own facilities for easy access. It is easy-to-use, requires only short hands-on time and provides quick and reliable results.

Serial testing solutions **tailored** to your organisation

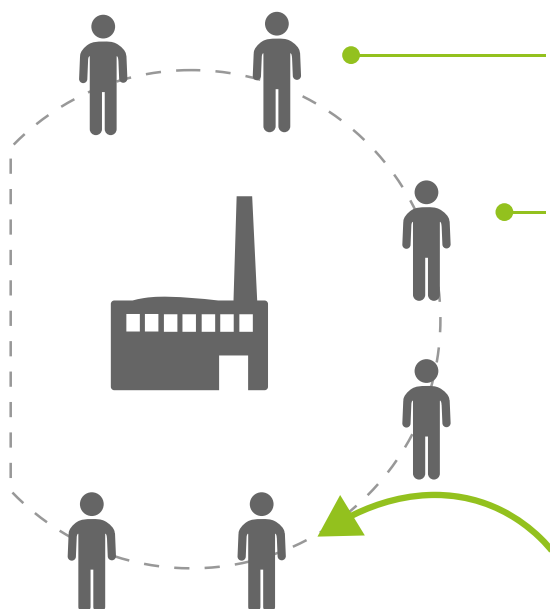

Test each employee few times per week

Suitable for factories, shipyards, construction sites, paper mills, oil rigs, etc.

Detect positive cases at the earliest

By reducing the time of circulation of a contagious individual, the number of contacts are lower, and less employees have to be isolated.



300 people tested in 24h

Even 300 workers can be tested in 24 hours with only 1 analyser. This way you can keep the company COVID-19 free.